

LATEST 360° CRYOLIPOLYSIS MACHINE

SEMICONDUCTOR COOLING

- More Comfortable
 Twice Treatment Areas
 Say Goodbye To Cellulite






MACHINE FUNCTIONS

Suitable for all skin types







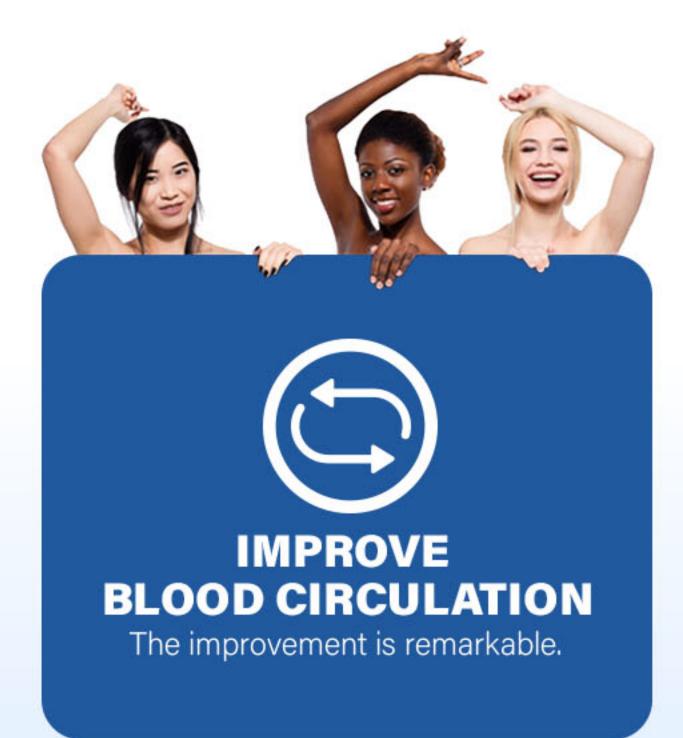



ADVANTAGES

COMFORTABLE

Medical level material, more smooth and soft.

HIGH EFFICIENCY

Multiple areas can be treated at one time.

SAFE

Vacuum increase gently, no harm to skin.

MORE INTELLIGENT

3 minutes heating massage, enhance blood circulation

EASY

Easily replace probe, Easy operation.

BETTER RESULTS

Full-scale freezing, larger treatment area.

IMMEDIATE EFFECT

Reduce size by 3cm in one treatment.





4 HANDLES 6 REPLACEMENT HEADS














6 HEADS

4 handles can works at same time and independently, 6 different size replacement heads can for different treatment areas, which is more convenient and saving times.



TREATMENT HEAD 1

TREATMENT AREA: 110×50.5MM

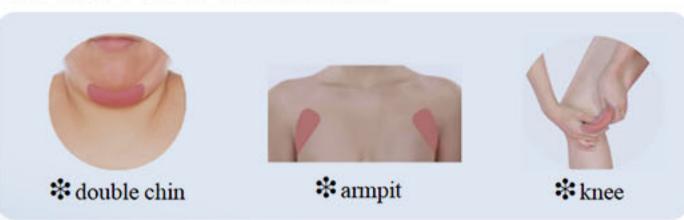


Suitable for small fat pockets and curved contours



TREATMENT HEAD 2

TREATMENT AREA: 110×80MM

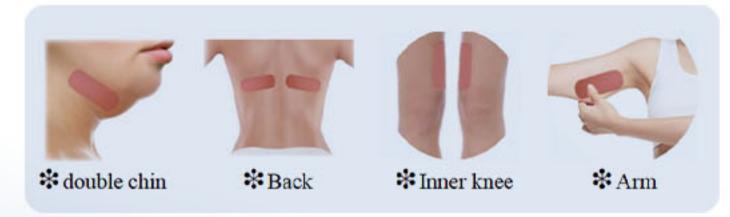


Suitable for small fat pockets and flat head contours



TREATMENT HEAD 3

TREATMENT AREA: 260×140MM

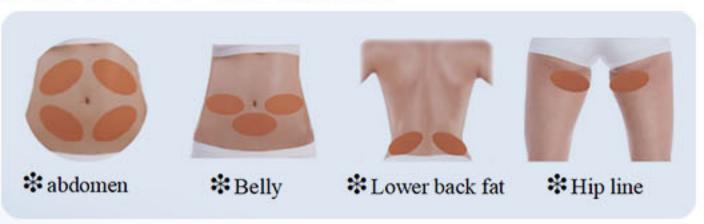

Suitable for wide and curved contours






TREATMENT HEAD 4

TREATMENT AREA: 330×180MM



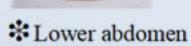



*Upper and lower abdomen



*Back fat

Suitable for large fat pockets and large fat deposits

TREATMENT HEAD 5


TREATMENT AREA: 200×140MM









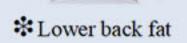




*Hip line

Suitable for large fat bags and flat contours



TREATMENT HEAD 6

TREATMENT AREA: 250×190MM

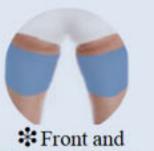


Front and back of thighs



*Outer thigh



* Waist fat

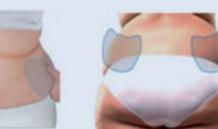

*Lower back

Suitable for the thigh area; flat and wide fat areas







360° GLAMOURESSENCE SURROUND COOLING




The 360° surround cooling technology can increase the efficiency by up to 18.1%, and the adsorbed area is frozen in a 360° omni-directional way without dead corners, resulting in more effective removal of fat cells.

360° ALL-AROUND COOLING:




Cooling energy complete and uniformity.

TRADITIONAL TWO-SIDE COOLING VERSION:

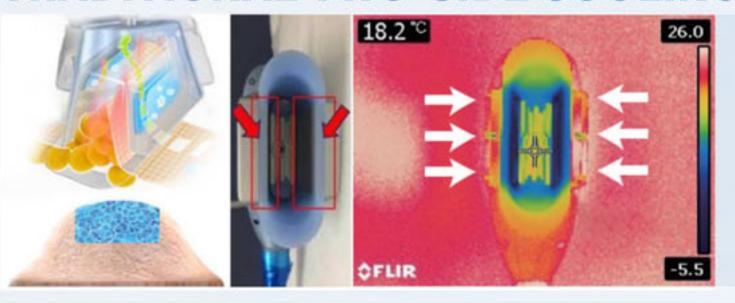


Cooling energy focus on both sides.

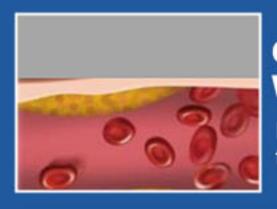

COOLING PLATE WITHOUT NEGATIVE

Blood flow and boady temperature tend to remain constant

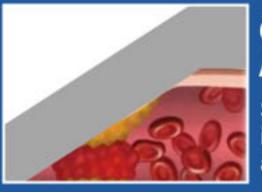

COOLING **AFTER ADSORPTION**

Suction blocks blood flow, allowing low temperature to reach adipose tissue.

The clinical effect is beetter if blood circulation is blocked and fat freezing is accelerate.






The general 360° cooling plate areain the market






TREATMENT AREA

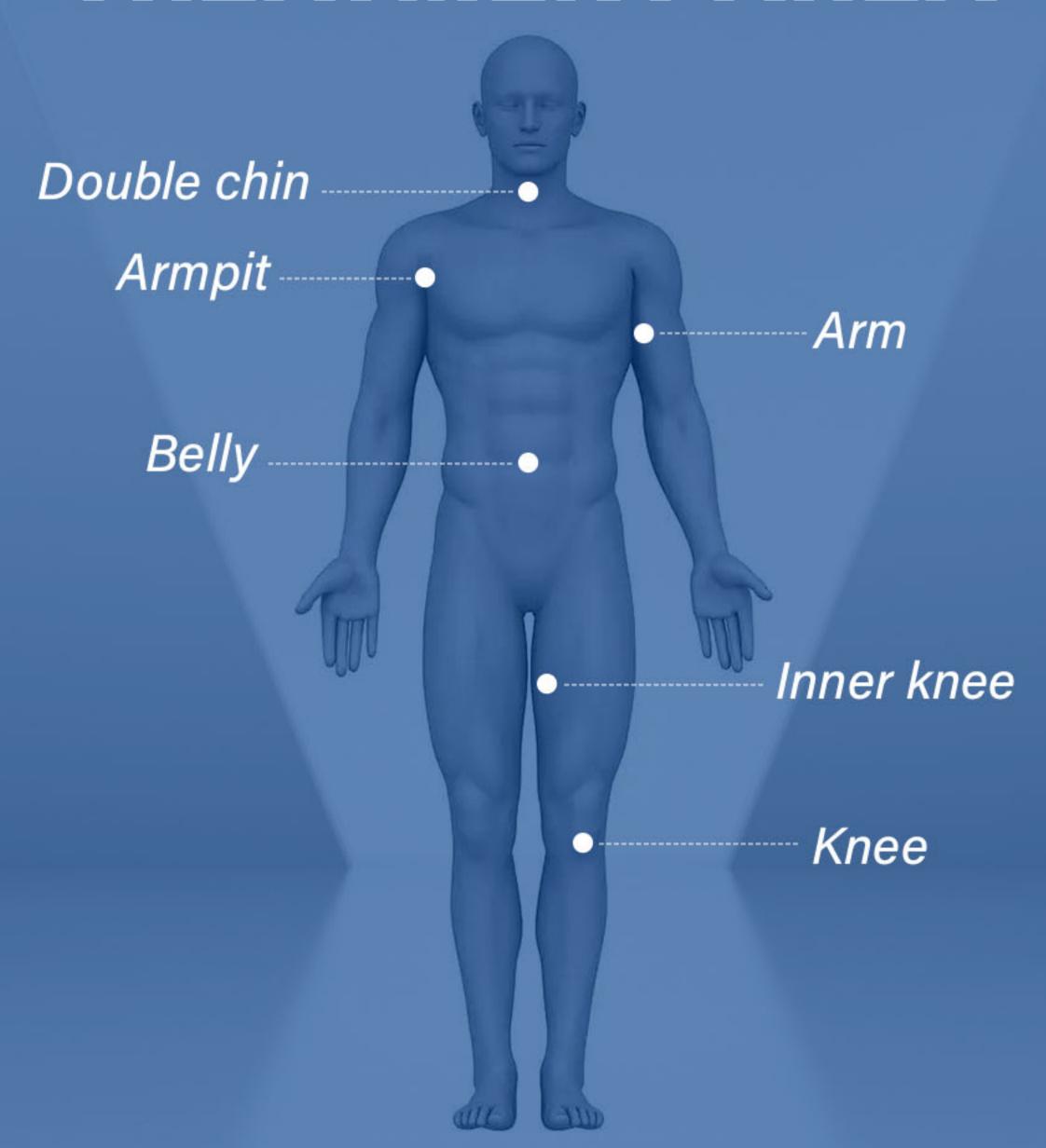



BEFORE AND AFTER

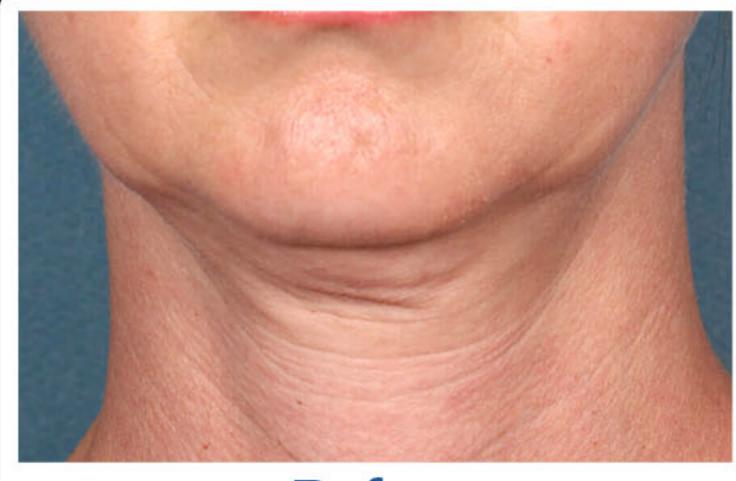


Before



After



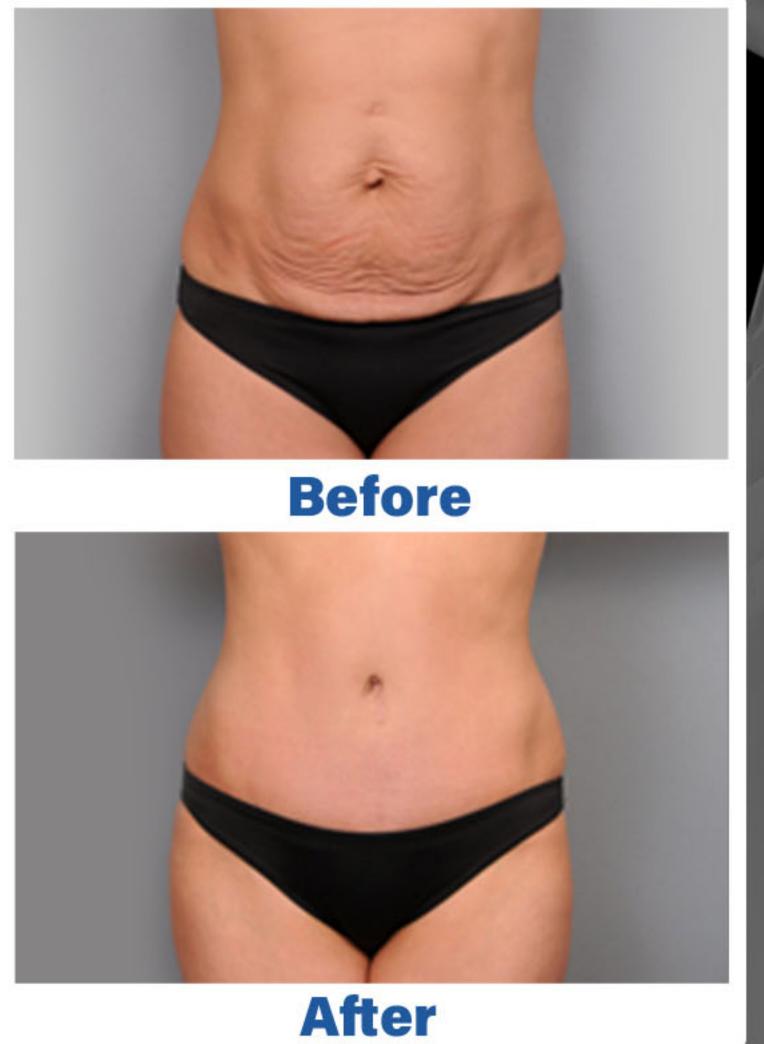






15.6 INCH **SMART SCREEN**



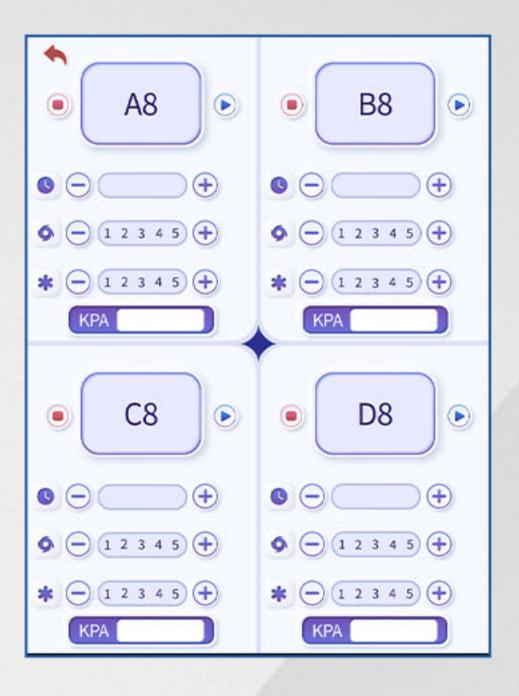

- Convenient operation Fast response time
- Encrypted design for privacy



MACHINE PARAMETERS

Machine Name Diamond Ice Body Slimming Machine

Display Screen 15.6 Inch Smart Screen

Cooling Temperature -15-50°C

Vacuum Suction 1-5 Gears (10-50Kpa)

Setting time 1-99min (Default 60min)

Output Power 1000W


Input voltage 110V/220V

Handle 4 Handles (6 Heads)






WORKING PRINCIPLE

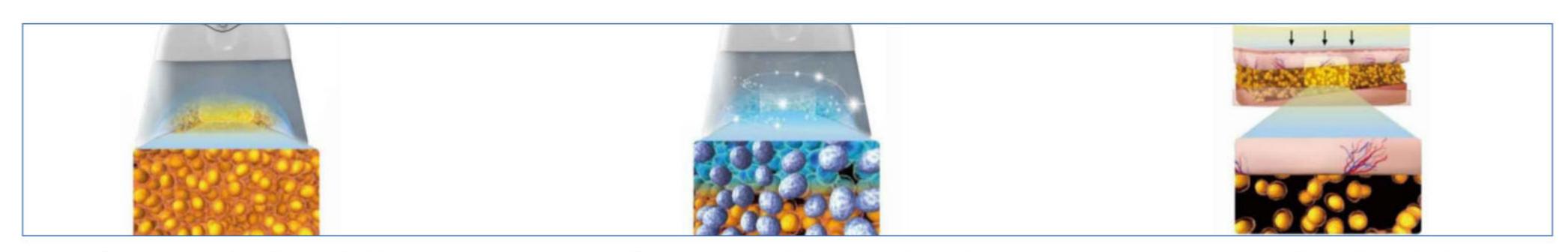

The treatement head is attached

The cooling energy treats the target fat layer without damaging

Damaged fat celles are broken down and cleared throught the natural body process(apoptosis)

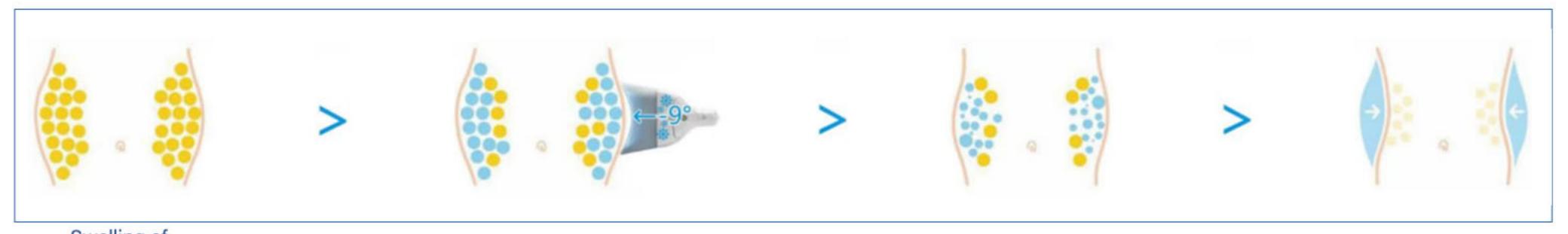

Swelling of abdominal fat

Convey cooling

fat cells are destroyed

decreased number of fat cells in the abdomen